

Otto Bock®

QUALITY FOR LIFE



Manu Forsa Volar

High stability for the wrist

The 50P12 Manu Forsa Volar wrist orthosis provides for a high degree of immobilization of the wrist in a functional position.

It is made out of a skin-friendly velour material and is fabricated using state-of-the-art fabrication technologies. A flexible plastic reinforcement insert and a volar, Y-formed, anatomically pre-shaped metal splint are integrated into the orthosis to stabilize and support the wrist. A heart-shaped pad, attached by hook and loop closure, supports the palm.

Three circular hook and loop straps allow for continuous individual adjustment, and make the wrist orthosis easy to apply. The orthosis' appealing design also contributes to patient compliance.

Otto Bock HealthCare GmbH

Max-Näder-Straße 15 · 37115 Duderstadt/Germany
Phone +49 5527 848-0 · Fax +49 5527 848-1414
healthcare@ottobock.de · www.ottobock.com



Manu Forsa Volar

50P12 Wrist Orthosis

Indications


- Acute and chronic irritations in the wrist (such as carpal tunnel syndrome)
- Post-operative or post-traumatic conditions, which require a high degree of immobilization of the wrist
- Partial paralyses

Effects

The 50P12 Manu Forsa Volar wrist orthosis provides for a high degree of immobilization of the wrist in a functional position.

Features


- The functional Y-splint provides stable, volar support.
- For comfortable, anatomical positioning, the palm rests on a pad which is attached by hook and loop closure.
- Skin-friendly, soft velour material for high comfort
- Circular hook and loop straps for easy adjustability and application
- Low weight
- Attractive design

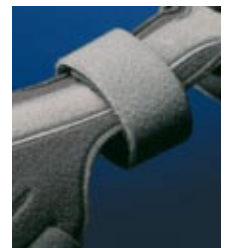

Size selection

The 50P12 Manu Forsa Volar wrist orthosis is available in two sizes (for right and left hands):

Order no.	Side	Size	Hand width (A)
50P12=R-S	right	S	up to 8.5 cm
50P12=L-S	left	S	up to 8.5 cm
50P12=R-M	right	M	from 8.5 cm
50P12=L-M	left	M	from 8.5 cm

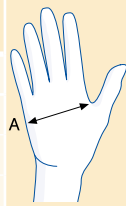

Sender

Order no.

Quantity

50P12=R-S

pc(s)

50P12=L-S

pc(s)

50P12=R-M

pc(s)

50P12=L-M

pc(s)



Customer ref.

Customer no.

Otto Bock[®]

QUALITY FOR LIFE