


## ATMOSAIR® 9000 MRS


The only non-powered, clinically and cost-effective pressure redistribution surface with SAT™



# ATMOSAIR® 9000 MRS

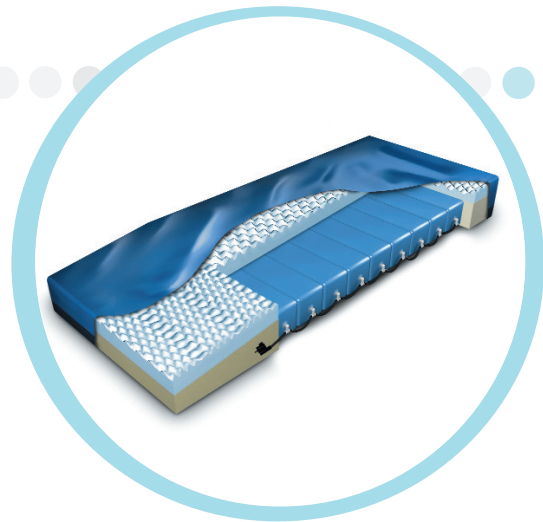

SAT™ system design

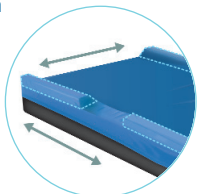

Feature	Benefit
Self Adjusting Technology™ (SAT™): A system of dynamic air chambers and cut-off valves	<ul style="list-style-type: none"> <li>Helps to maximize body weight displacement and minimize tissue interface pressure by automatically adjusting the internal air pressure in reaction to body movement, altitude and temperature.</li> </ul>
9 Air Cylinders	<ul style="list-style-type: none"> <li>Support under the shoulders, seat and legs via nine zones of horizontal pressure distribution.</li> </ul>
Non-Powered Pressure Redistribution	<ul style="list-style-type: none"> <li>Uninterrupted pressure relief even in the event of an electrical power failure.</li> </ul>
Sloped Heel Section	<ul style="list-style-type: none"> <li>Helps to transfer weight off of the heel and onto the calf and thigh area of the leg.</li> </ul>
Firm Perimeter	<ul style="list-style-type: none"> <li>Provides seating support for the caregiver or patient.</li> </ul>
Raised Side Bolsters (optional)	<ul style="list-style-type: none"> <li>Serve as a passive reminder of the mattress edges and facilitate patient ingress/egress.</li> </ul>
Customizable Length and Width (optional)	<ul style="list-style-type: none"> <li>Helps ensure proper fit on a wide variety of bed frames.</li> </ul>

## Quality Materials

### Cover

- TheraTex™ Nylon: A durable fabric with a blended/butyl rubber backing and a low MVT rate. (For use primarily in higher turnover Critical Care areas.)

### Bottom

- SoffTICK® Low Slip: A 10 oz. non-skid bottom material. (For use primarily on smooth, pan-type bed frame.)
- Centurion®12: A 12 oz. vinyl bottom material. (For use primarily on spring-type bed frame.)

## Ordering

KCI Order Number	Description
KA-9-SONB-SIZE**	AtmosAir® 9000, TheraTex® Top, Non-skid Bottom
KA-9-SOVG-SIZE**	AtmosAir® 9000, TheraTex® Top, Vinyl Gray Bottom

\*\*Please refer to Product Specifications for sizing, or contact your local KCI sales consultant regarding custom sizing options.

Individual facility case studies have documented reduced incidence of pressure ulcers and reduced costs through purchase and use of AtmosAir® mattresses as part of a wound care program.\*\*

For more information, visit [www.kci-medical.com](http://www.kci-medical.com)

**KCI THERAPEUTIC**  
Support Systems

### Product Specifications\*

Length . . . . .75"-84" by 1"  
 Width . . . . .32"-37" by 0.5"  
 Depth . . . . .7"  
 Optional Bolsters . . . . .2" high with 18" cutout  
 Anti-microbial: Inherent in both the cover and foam.  
 Flammability: All mattress materials meet TB 117 standards and all mattresses meet Federal cigarette ignition standards.  
 TB 603, ASTM E 1590, TB 129 or BFD IX-11  
 Compliance: Available with an optional fire sleeve at additional cost.

Mattress Warranty . . . . .5 years non-prorated  
 Maximum Patient Weight . . . . .500 lbs.†  
 Air Pump (optional)  
 Warranty . . . . .2 years non-prorated  
 Power . . . . .110V / 60Hz  
 Output . . . . .8 liters per minute  
 Cycle Time . . . . .10 minutes  
 Weight . . . . .5 lbs.

\*Specifications made to order. Allow 6 weeks for delivery.  
 †Consult frame manufacturer regarding individual frame capacity.

\*\*Case studies are based on single facility findings. Individual results may vary. Studies available upon request.

AtmosAir® MRS units have specific indications, contraindications, safety information and instructions for use. Please consult product labeling and instructions for use prior to use.

©2008, KCI Licensing, Inc. All rights reserved. KCI USA, Inc., PO Box 659508, San Antonio, Texas 78265-9508, 1-800-275-4524. In Canada, call 1-800-668-5403. SoffTICK® is a trademark of Vintex.® All other trademarks and service marks designated herein are property of KCI, its affiliates and licensors. Most KCI products are patented and/or subject to pending patents. KCI is an ISO-13485:2003 certified company. Rev 08/08 Lit.# 61-A-008

